



The design concept is preventing troubles with a simple structure.

Tierrett Ureter Bypass Tube (Soft type)

Registration No. of Medical Equipment for Animal Use: 1 Douyaku 447



This product is placed in the kidney and bladder to assist the urinary function of the urinary tract.

Specifications

Size		Tip length	Effective length	O.D.	I.D.	Tip side hole	Product cord	JAN cord
8.5 Fr	For cats	35 mm	500 mm	2.8 mm	1.9 mm	3 holes	800-001-7154	4542187228205
	For dogs	40 mm					800-001-7155	4542187228212
10 Fr	For cats	35 mm		3.3 mm	2.1 mm		800-001-7152	4542187228182
	For dogs	40 mm					800-001-7153	4542187228199

1-pack box / EOG sterilized

Please note that the product specifications are subject to change without notic to reflect the results of our consist R&D efforts.

Read the package insert before using the product.




Product features

(5 ideas to keep the urinary tract safe)

1 The renal pelvis insertion part is pre-bent in a L-shape

to prevent obstruction of the lumen due to bending.

- 2 The loop shape absorbs body movements, reduces the burden on the renal pelvis, sutured part of the bladder, etc., and prevents a tube coming off.
- 3 A thick tube (8.5 Fr: 1.9 mm of I.D., 10 Fr: 2.1 mm of I.D.) prevents clogging.
- 4 A tube has no joints so urine flows without stagnation.
- ⑤ No accesses from outside the body provids a low risk of infection.

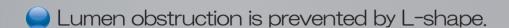

Tip length For cats: 35 mm For dogs: 40 mm Absorbing body movements leads to reduce the burden.

2 pcs of fixing plates


(They are attached to a tube. Please use medical instant adhesive for fixing.)



Protective sheets

Shape protection wire (It keeps a tube straight.

Pull out and discard.)

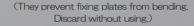



CREATE MEDIC CO., LTD.

https://www.createmedic.co.jp/english/

(Contact information)

Animal Medical Business Department E-mail: oem-infoch@createmedic.co.jp

[Head office] 2-5-25 Chigasaki Minami, Tsuzuki-ku, Yokohama, Kanagawa, Japan



Control No: 01-0447-20210624